



RESTORE
PEACE-OF-MIND
AND ACCESS
IN YOUR HOME

www.bruno.com



COMBINE QUALITY, COMFORT AND STYLE

When you want the highest quality stairlift combined with a stylized appearance, you want the Bruno Elite. The Elite offers worry-free access to all levels of your home. And the Elite's stylized design blends with your home's environment.

Let Bruno's Elite safely connect you to all levels of your home again.

- Made-in-the-USA quality
- Easy, reliable operation
- Smooth, quiet ride
- Luxury styling and comfort

400 lb Lift Capacity (190kg)

Learn more at www.bruno.com.

BRUNO HAS AN INDUSTRY-LEADING WARRANTY

Gain peace of mind with the most comprehensive manufacturer warranty in the stairlift industry.

Lifetime Warranty on Major Components -- motor, gearbox, drivetrain and rail -- for the lifetime of the Bruno Stairlift as long as the original purchaser owns the stairlift.

2-Year Warranty on All Components after date of installation by an Authorized Bruno Dealer.

30-Days Labor Coverage after date of installation by an Authorized Bruno Dealer for any manufacturing defect.

Ask your Bruno dealer for full details!





ELITE FEATURES

- Plush, generous-size seat with multiple height adjustments
- Adjustable width armrests
- Anodized rail with covered gear rack gives a sleek appearance
- Flip up arms, seat, and footrest create plenty of extra space on steps
- Adjustable footrest for maximum comfort
- Vertical rail design enables tight installation to within 6 in (152 mm) of wall
- Onboard audio/visual diagnostics for easy service
- Continuous charge strips along rail recharge unit wherever it stops
- Two wireless call/sends allow unit to be controlled remotely
- Two 12v batteries ensure dependable performance even in a power outage
- Soft start/stop gives a smooth ride from start to finish




Ask your local dealer about upholstery and color options.



Offset swivel seat makes it easy to get on/off.

“The Bruno stairlift has eliminated an enormous amount of anxiety and stress. We feel grateful for it every day.”

A. Wambold - VA

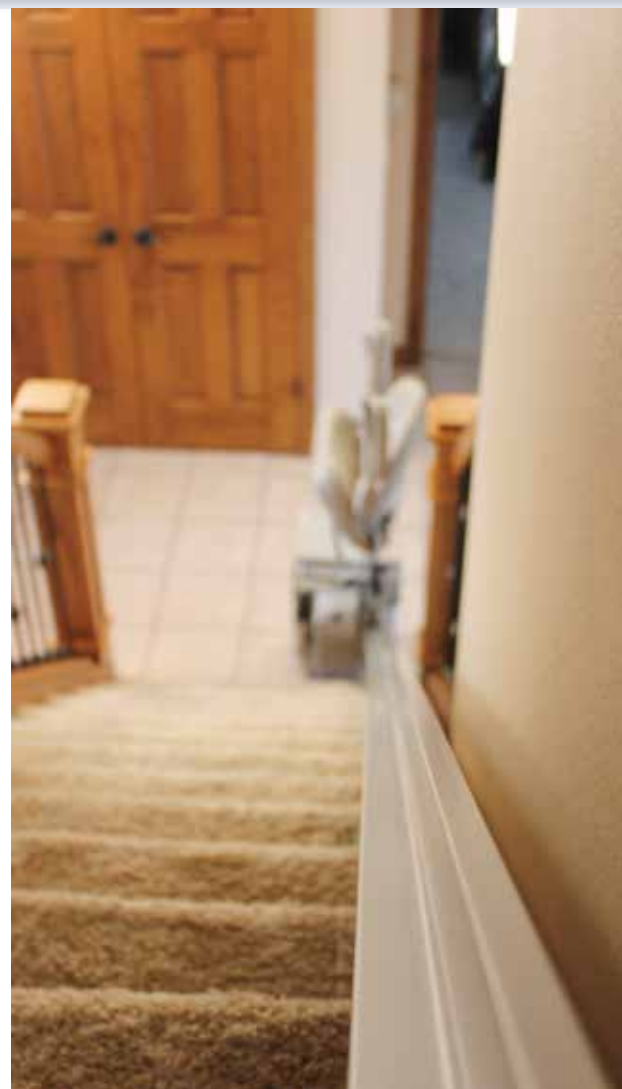


Close installation to wall.

ELITE OPTIONS



Power folding footrest automatically flips up/down when seat is raised/lowered.

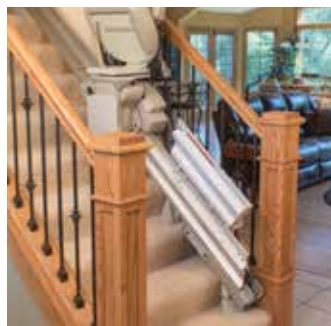

Power or manual folding rails for narrow hallway or when doorway is at the bottom of stairs. Manual operation or push button automatic.



Power swivel seat for effortless exit. Controlled on chair arm or wireless call/send.



Larger seat and footrest for increased comfort.



"The stair lift is great, very solid and looks good too."

L. UHRING - MI



TECHNICAL SPECIFICATIONS

Standard Seat Size: 17.5 in wide x 15 in deep
(445 mm wide x 381 mm deep)

Larger Seat Option: 21.5 in wide x 16 in deep
(546 mm wide x 406 mm deep)

Standard Footrest Size: 10 in long x 9 in wide
(254 mm long x 229 mm wide)

Larger Footrest Option: 12 in long x 11 in wide
(305 mm long x 279 mm wide)

About Bruno

Bruno Independent Living Aids has helped improve the lives of people with limited mobility for nearly over 35 years. Engineered and manufactured in the USA, Bruno stairlifts and vertical platform lifts help people better access homes, public buildings and businesses, and Bruno's full line of scooter/powerchair lifts help people easily transport their mobility devices. A family-owned business, Bruno takes pride in their unwavering focus on quality and safety and their position as industry innovation leader.

Proudly made in the USA.



Bruno Independent Living Aids, Inc.
1780 Executive Drive, Oconomowoc, WI 53066

www.bruno.com

© Bruno Independent Living Aids, Inc. 2019
2018-03-20-Elite-6-PageV4.indd • MKT-SL-2



CONSUMERAFFAIRS